

Serum Collection Set

The device is used for preparing autologous serum in a completely closed system by simple operation

FEATURE

◆ Closed system

Serum can be prepared, separated, and stored in a completely closed system. There is no risk of bacterial contamination.

◆ Simple operation

Simple operation: blood collection, agitation, and centrifuging. Enable to prepare serum in approximately 1 hour.


◆ Platelet activation beads

When platelet activation beads contact with the blood drawn into the collection bag, they promote blood coagulation, thus enabling the separation of serum in a relatively short time.

SPECIFICATION

Item	Bag Type	Tube Type
Product code	902230000	909530400
Volume of blood drawn	70~200 mL	15~20 mL
Sterilization method	Steam sterilization	
Expiry Date	3 years	
Details	<ul style="list-style-type: none"> • Collection bag – 200 mL • Spare bag – 150 mL • Storage bag – 50 mL × 3 	<ul style="list-style-type: none"> • Collection container – 20 mL • Storage container – 10 mL
Package	1 Set × 5 Boxes = 5 Sets	4 Sets × 2 Boxes = 8 Sets

Large volume type for convenient preparation of a large amount of autologous serum

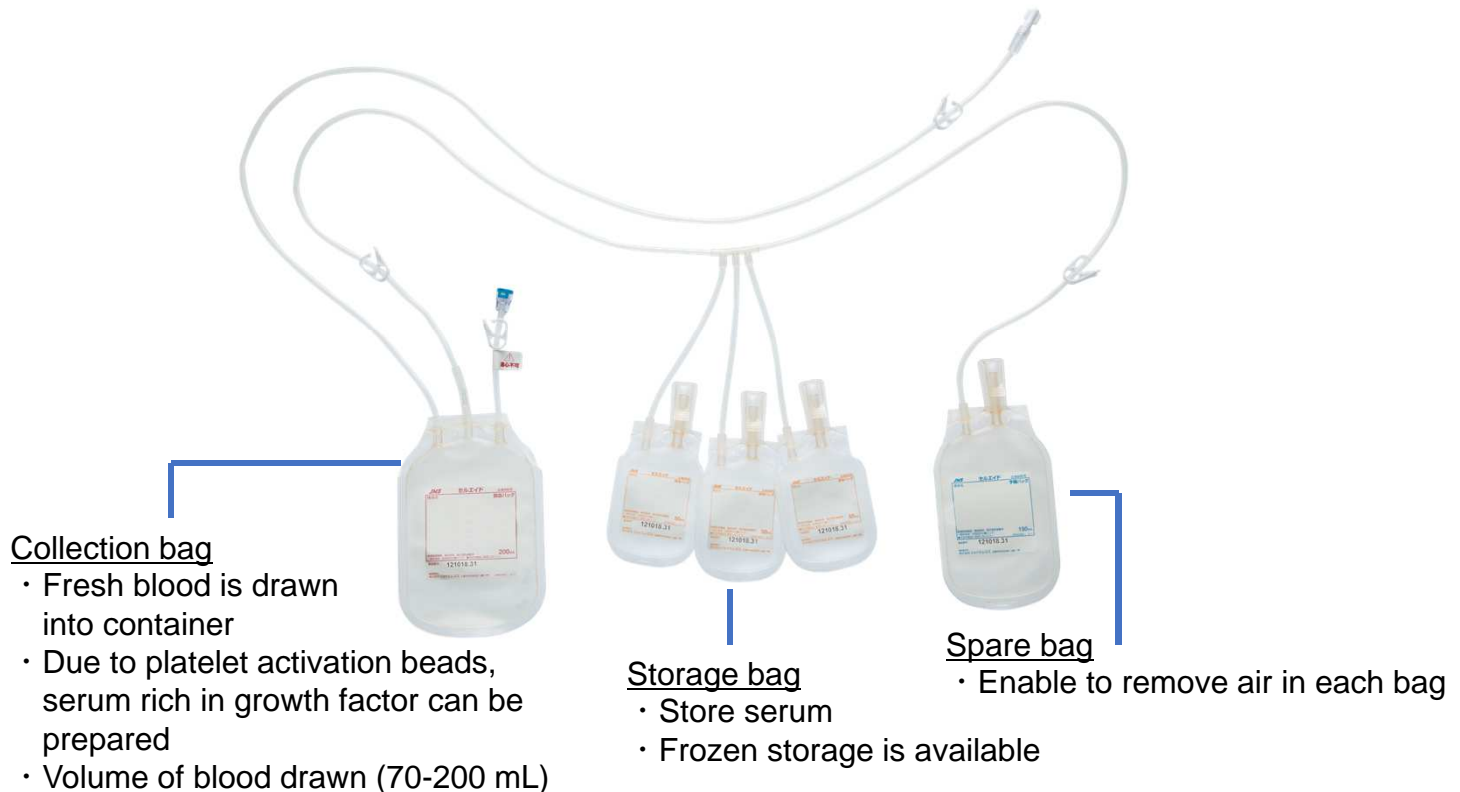

Compact type for small volume preparation of autologous serum

