



A.V. FISTULA NEEDLE SET

High tech with decades of experience




TWIN PACK

JMS recommends the use of a back eye needle for the arterial line and a without back eye for the venous line. JMS's Twin Pack (one arterial and one venous needles in a blister pack) offers stocking convenience and economy. Color coded red and blue clamps specify the lines and simply eliminate the risk of "reversed line" recirculation.



SINGLE NEEDLE SET

This set includes a carefully designed Y - piece joint and rotating hub. This combination permits the convenience of needle bevel rotation and secure compression free fixation as well as high dialysis efficiency by using minimal dead space which allows for a smooth flow without turbulence.

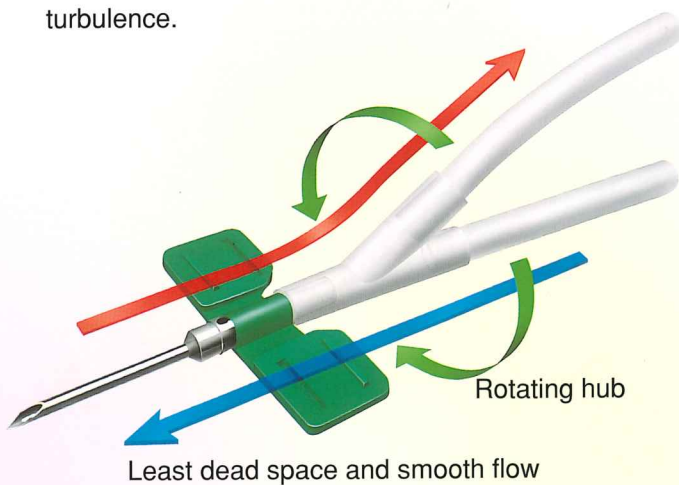

Packaging features :

- Peel-to open blister package for proper sterile handling.
- Color coded and symbolized labeling on the shelf box.
- Choice of convenient dispenser opening on the side panel or fully opening top panel.



Other applications

JMS fistula needles are used mostly for arterio-venous fistula or graft in hemodialysis treatment, but are also widely used as blood access devices for other extracorporeal applications. Its laminar flow inner configuration without gap and siliconized inner lumen greatly reduce clot formation and damage to components. JMS needles are much appreciated for apheresis and automated blood components collection, where blood flow is low and higher quality of blood component are demanded.

JMS



Others

